

The family of One Solution st™ Water For Injection products meets United States Pharmacopeia (USP) Water For Injection (WFI) specifications and is intended to be used to dilute surface disinfectants and cleaners. Sterile filtered and filled in a controlled environment, each package is double bagged and gamma irradiated for easy entry into your controlled environments and aseptic core areas. These products have an extensive documentation package to meet your validation objectives and are supported by the STERIS Technical Services team for individualized application and validation assistance.



## One Solution st™ Water For Injection

High quality water to dilute surface disinfectants and cleaners

### Features

Meets USP Water For Injection specifications

Produced by Finn-Aqua® Multiple Effect Water Still

Double bagged, 0.2 micron filtered, gamma irradiated and each lot USP sterility tested

Manufactured under a FDA cGMP, ISO 9001 and ISO 13485 compliant quality system

### Benefits

Extremely low level of endotoxins reduces contamination concerns

Highest quality Water For Injection for your facility


Allows easy entry into controlled environments

Supports quality assurance and quality control

### Applications

Dilution of concentrated cleaner or disinfectants

- Residue removal from:
  - » Floors
  - » Walls
  - » Ceilings
  - » Countertops
  - » BSL Hoods
  - » Carts and Other Equipment
- Products may be applied by:
  - » Mopping
  - » Wiping
  - » Spraying

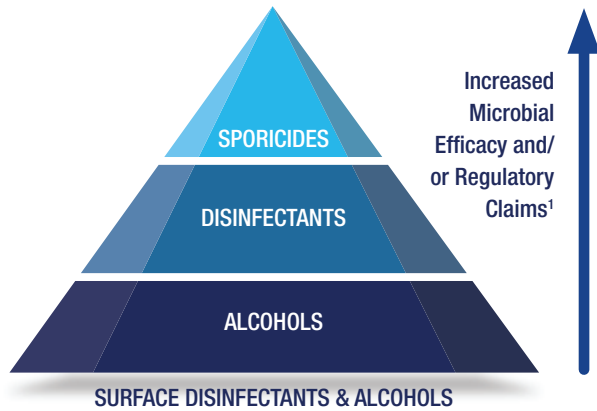

<sup>1</sup> Products that fall into the categories at the bottom of the pyramid are most frequently used and are generally not sporicidal. Progression up the pyramid indicates stronger performance overall and a broader spectrum of claims.

In maintaining microbial control of your facility, what is used to dilute disinfectants can be just as important as the disinfectants themselves. Regulations require the reduction of contaminants from all solutions introduced into your controlled environments or aseptic core areas. Similarly, when removing residues from product contact surfaces and the rest of your controlled environments, it is important to limit introduction of potential contamination and residues from the cleaning solution. One Solution st Water For Injection is ideal for these important applications because it is manufactured to minimize contaminants with no significant residues. In addition, our Technical Service group can help customize your contamination control program to take advantage of this excellent product.

### One Solution st™ Water For Injection

High quality water to dilute surface disinfectants and cleaners

### Ordering Information (Europe)

Description	Part Number
1 gallon (3.78 L) bottle (4 per case)	1S0308WR
2 gallon (7.57 L) bottle (2 per case)	1S03E4WR