

Septihol Ethyl WFI Sterile Alcohol Solution Products are a mixture of 70% (v/v) reagent-grade alcohol and 30% (v/v) United States Pharmacopeia (USP) Water For Injection (WFI). They are designed to remove organic and inorganic residue from product and non-product contact surfaces. Sterile filtered and filled in a controlled environment, each package is double bagged for easy entry into your controlled environments and aseptic core areas. These products have an extensive documentation package to meet your validation objectives and are supported by the STERIS Technical Services team for individualized application and validation assistance.



Septihol™ Ethyl WFI Sterile Alcohol Solution

Ethanol based cleaner

Features

Diluted with USP Water For Injection

Formulation easily evaporates

Double bagged

0.2 micron filtered, filled in a controlled environment, gamma irradiated and each lot USP <71> sterility tested

Manufactured under a FDA cGMP, ISO 9001 and ISO 13485 compliant quality system

Benefits

Supports regulatory compliance and endotoxin control goals

Appropriate for use on your most critical surfaces


Allows easy entry into critical environments

Facilitates introduction into critical environments

Supports quality assurance and quality control goals

Applications

- Residue removal from:
 - » Product contact Surfaces
 - » Floors
 - » Walls
 - » Ceilings
 - » Countertops
 - » BSL Hoods
 - » Carts and Other Equipment
 - » Gloves
- Products may be applied by:
 - » Mopping
 - » Wiping
 - » Spraying

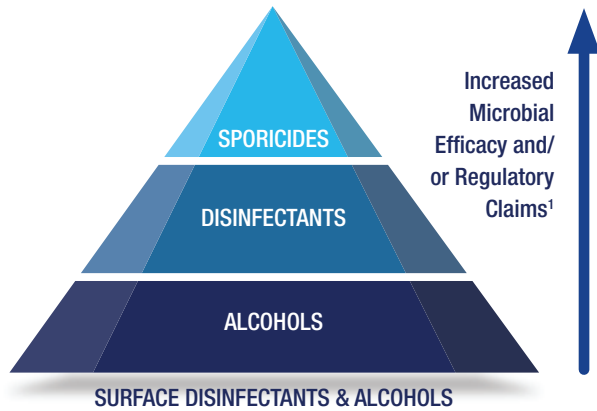

¹ Products that fall into the categories at the bottom of the pyramid are most frequently used and are generally not sporicidal. Progression up the pyramid indicates stronger performance overall and a broader spectrum of claims.

An important aspect to any microbial control program is including rinsing or cleaning steps to control residues. When removing residues, it is important to use products that do not leave behind a different residue. The only ingredients in our Septihol Ethyl WFI Sterile Alcohol Solutions are ethanol, isopropanol, methanol and WFI. This makes Septihol Ethyl WFI Sterile Alcohol Solutions ideal for this important application since it easily evaporates. Our Technical Services team can help customize your contamination control program to take advantage of this excellent product.

Septihol® Ethyl WFI Sterile Alcohol Solution

Low residue ethanol based cleaner

Ordering Information (Europe)

Description	Part Number
22 oz (650 mL) trigger sprayer (12 per case)	1S1077WR
32 oz (946 mL) trigger sprayer (12 per case)	1S1096WR
1 gallon (3.78 L) bottle (4 per case)	1S1008WR